

# TNFAIP8L2 Antibody


## Basic information:

<b>Catalog No.:</b>	UPA03849	<b>Source:</b>	Rabbit
<b>Size:</b>	50ul/100ul	<b>Clonality:</b>	Polyclonal
<b>Concentration:</b>	1mg/ml	<b>Isotype:</b>	IgG
<b>Purification:</b>	The antibody was purified by immunogen affinity chromatography.		

## Useful Information:

<b>Applications:</b>	WB (1:500 - 1:1000)
<b>Reactivity:</b>	Human, Mouse, Rat, Pig, Rabbit
<b>Specificity:</b>	Recognizes endogenous levels of TNFAIP8L2 protein.
<b>Immunogen:</b>	KLH-conjugated synthetic peptide encompassing a sequence within the center region of human TNFAIP8L2. The exact sequence is proprietary.
<b>Description:</b>	Rabbit polyclonal antibody to TNFAIP8L2
<b>Uniprot:</b>	Q6P589
<b>BiowMW:</b>	Refer to Figures
<b>Buffer:</b>	Liquid in 0.42% Potassium phosphate, 0.87% Sodium chloride, pH 7.3, 30% glycerol, and 0.01% Thimerosal.
<b>Storage:</b>	Shipped at 4°C. Upon delivery aliquot and store at -20°C for one year. Avoid freeze/thaw cycles.
<b>Note:</b>	For research use only, not for use in diagnostic procedure.

## Data:

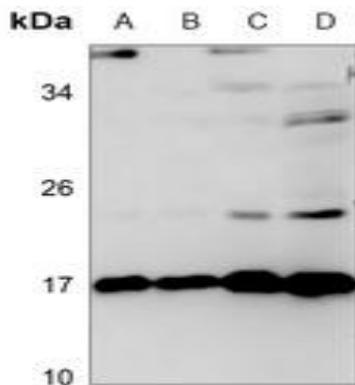

Western blot analysis of TNFAIP8L2 expression in mouse liver (A), mouse kidney (B), rat liver (C), rat kidney (D) whole cell lysates. (Predicted band size: 20 kD; Observed band size: 17; 21 kD)