

p53 (Phospho-T18) Antibody


Basic information:

Catalog No.:	UPA03854	Source:	Rabbit
Size:	50ul/100ul	Clonality:	Polyclonal
Concentration:	1mg/ml	Isotype:	IgG
Purification:	The antibody was purified by immunogen affinity chromatography.		

Useful Information:

Applications:	WB (1:500 - 1:1000)
Reactivity:	Human, Monkey, Pig, Rabbit
Specificity:	Recognizes endogenous levels of p53 protein only when phosphorylated at T18.
Immunogen:	KLH-conjugated synthetic phosphopeptide corresponding to residues surrounding T18 of human p53 protein. The exact sequence is proprietary.
Description:	Rabbit polyclonal antibody to p53 (Phospho-T18)
Uniprot:	P04637
BiowMW:	Refer to Figures
Buffer:	Liquid in 0.42% Potassium phosphate, 0.87% Sodium chloride, pH 7.3, 30% glycerol, and 0.01% Thimerosal.
Storage:	Shipped at 4°C. Upon delivery aliquot and store at -20°C for one year. Avoid freeze/thaw cycles.
Note:	For research use only, not for use in diagnostic procedure.

Data:

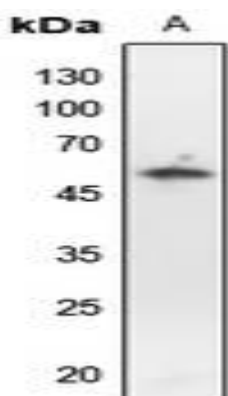

Western blot analysis of p53 (Phospho-T18) expression in HEK293T (A) whole cell lysates. (Predicted band size: 43 kD; Observed band size: 53 kD)