

TCP1 eta Antibody


Basic information:

Catalog No.:	UPA07315	Source:	Rabbit
Size:	50ul/100ul	Clonality:	Polyclonal
Concentration:	1mg/ml	Isotype:	IgG
Purification:	The antibody was purified by immunogen affinity chromatography.		

Useful Information:

Applications:	WB (1:500 - 1:2000)
Reactivity:	Human, Mouse, Rat
Specificity:	Recognizes endogenous levels of TCP1 eta protein.
Immunogen:	Recombinant full length protein of human TCP1 eta
Description:	Rabbit polyclonal antibody to TCP1 eta
Uniprot:	Q99832
BiowMW:	Refer to Figures
Buffer:	Liquid in 0.42% Potassium phosphate, 0.87% Sodium chloride, pH 7.3, 30% glycerol, and 0.01% Thimerosal.
Storage:	Shipped at 4°C. Upon delivery aliquot and store at -20°C for one year. Avoid freeze/thaw cycles.
Note:	For research use only, not for use in diagnostic procedure.

Data:

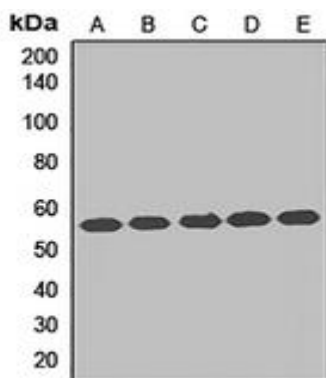

Western blot analysis of TCP1 eta expression in A431 (A), HeLa (B), NIH3T3 (C), mouse testis (D), rat testis (E) whole cell lysates. (Predicted band size: 37; 50; 54; 59 kD; Observed band size: 55 kD)