

ACADL Antibody


Basic information:

Catalog No.:	UPA07377	Source:	Rabbit
Size:	50ul/100ul	Clonality:	Polyclonal
Concentration:	1mg/ml	Isotype:	IgG
Purification:	The antibody was purified by immunogen affinity chromatography.		


Useful Information:

Applications:	WB (1:500 - 1:1000), IHC (1:50 - 1:200)
Reactivity:	Human, Mouse, Rat
Specificity:	Recognizes endogenous levels of ACADL protein.
Immunogen:	Recombinant full length protein of human ACADL
Description:	Rabbit polyclonal antibody to ACADL
Uniprot:	P28330
BiowMW:	Refer to Figures
Buffer:	Liquid in 0.42% Potassium phosphate, 0.87% Sodium chloride, pH 7.3, 30% glycerol, and 0.01% Thimerosal.
Storage:	Shipped at 4°C. Upon delivery aliquot and store at -20°C for one year. Avoid freeze/thaw cycles.
Note:	For research use only, not for use in diagnostic procedure.

Data:

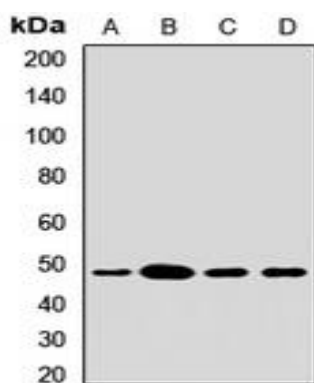

Western blot analysis of ACADL expression in SHSY5Y (A), mouse liver (B), mouse heart (C), rat spinal cord (D) whole cell lysates. (Predicted band size: 47 kD; Observed band size: 47 kD)

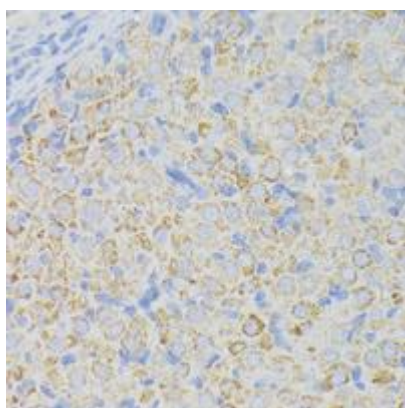

Immunohistochemical analysis of ACADL staining in rat ovary formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

